
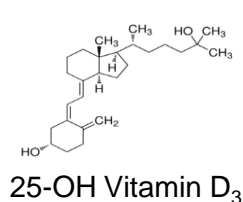
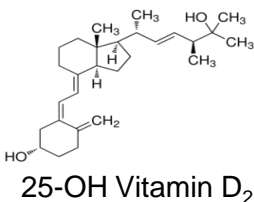


VITAMIN D METABOLITE EXTRACTION KIT


Compounds :

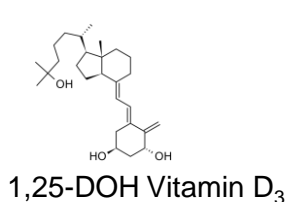

VITAMIN D KIT EXTRACTION METHOD:

(Vitamin D Metabolite Extraction kit: VDKT2C-30)

Automated Process using Oroflex Personal Pipettor

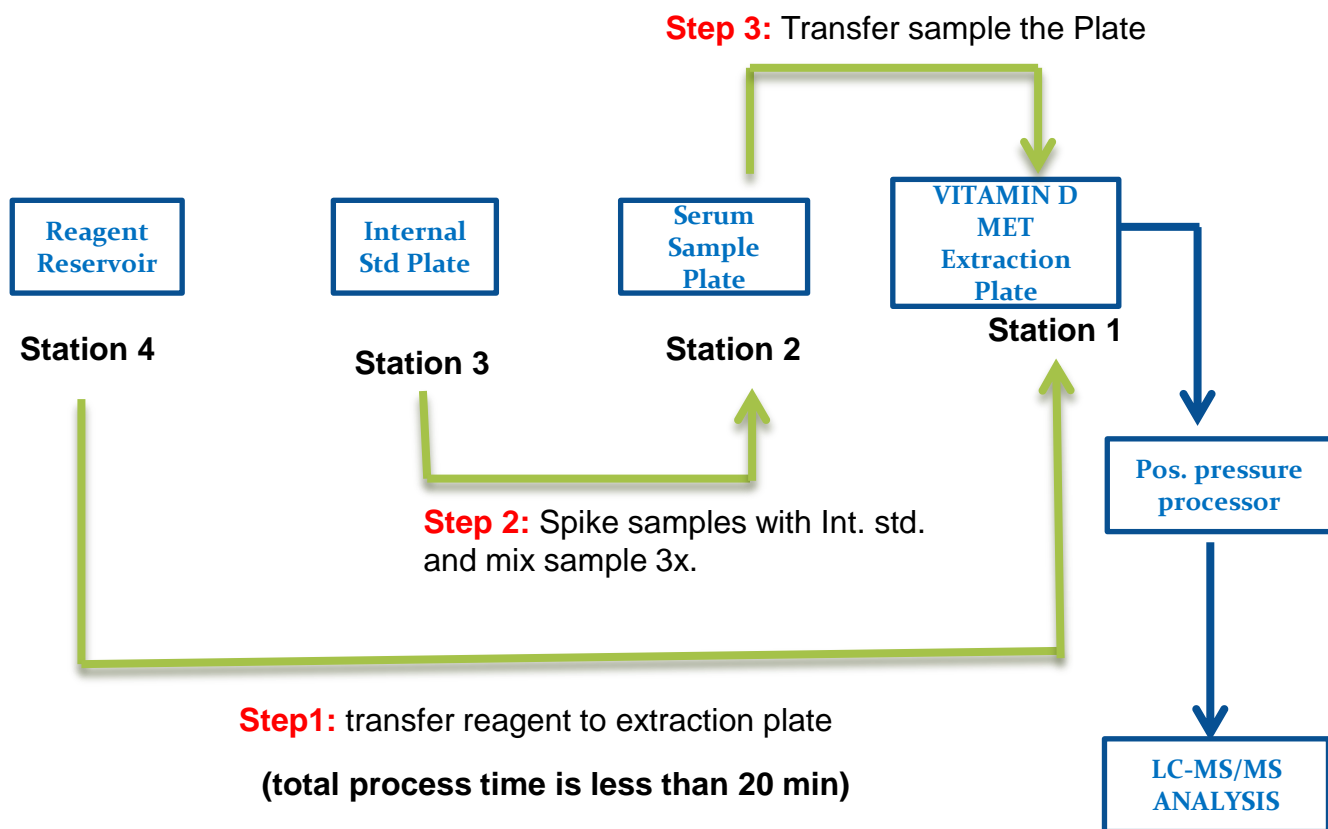

RESULTS:

	1,25-dihydroxy VD ₃	25-OH VD ₂	25-OH VD ₃
Recovery%	68	69	77
CV%	14	8	6

Phospholipid Profile for Serum Samples using Different Sample Preparation Products

Phospholipid profile

