


Vitamin D Metabolites Separated on Gazelle C18

Use: regulate calcium and phosphorus concentration in blood

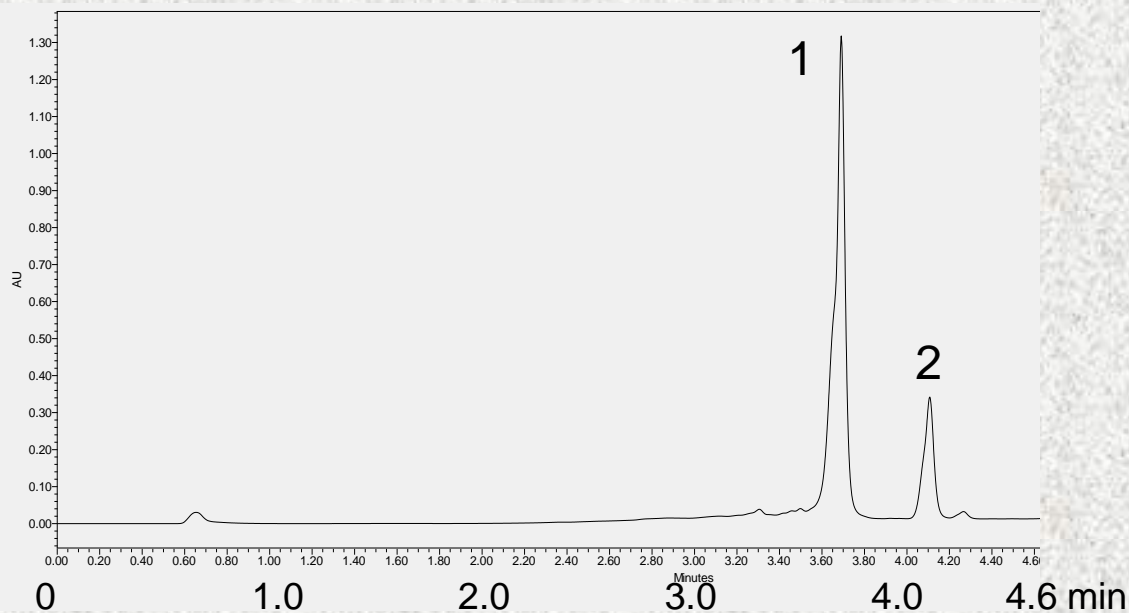

1	1,25-(OH) ₂ Vit D	2	25-OH Vit D
---	------------------------------	---	-------------

HPLC	Waters UPLC, H Class
Column	Gazelle C18, 1.7 μm
Column Size	2.1 x 50 mm
Temperature	25 °C
Detection	UV at 254 nm
Injection	10 μL

Mobile Phase	A: 15 mM NH ₄ -formate B: 15 mM NH ₄ -formate /MEOH		
Flow Rate	0.2 mL/min		
Gradient	Time	A	B
	0.00	60	40
	1.00	0	100
	4.60	0	100