



# Gazelle™

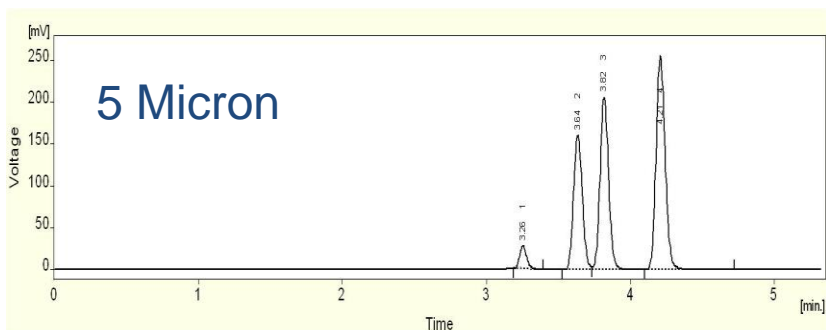

## 1.7 Micron UHPLC Amino Column

Orochem Technology's UHPLC columns offer the chromatographer significantly increased throughput, low solvent consumption, inline process monitoring, and improved LC-MS results.

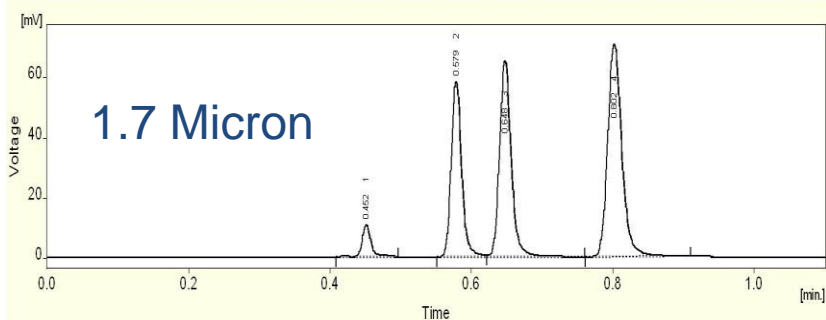

### Analytes

1. Naphthalene
2. Dipropyl phthalate
3. Diethyl phthalate
4. Dimethyl phthalate

	1.7 Micron	5 Micron
Column	3.0 x 50 mm	4.6 x 250 mm
Mobile Phase	97/3 Hept/IPA	97/3 Hept/IPA
Flow Rate	0.7 mL/min	1.0 mL/min
TP/m	166,600	80,200
Pressure	1650 psi	420 psi
Run Time	0.8 min	4.3 min