


1-Acenaphthenol on Epitomize CSP-1Z, 3 μ m, SFC-UHPLC

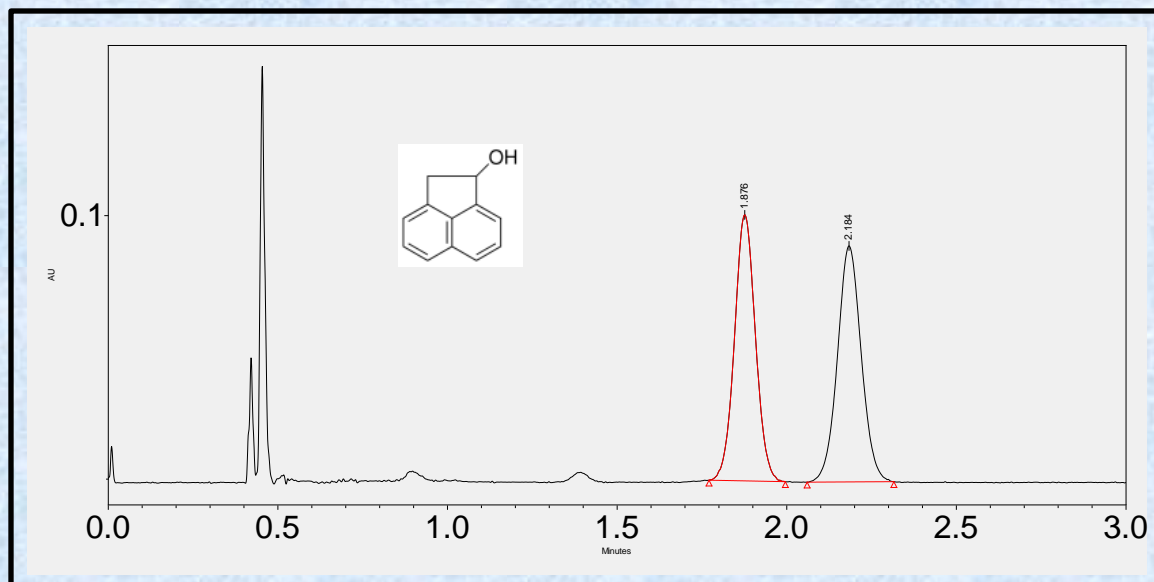

System	Waters Acquity UPC ²
Column	CSP-1Z, 3 micron
Column Size	3.0 x 50 mm
Mobile Phase	90/10 CO ₂ /Modifier
Modifier	50/50/0.2 MeOH/IPA/TFA
Flow Rate	0.80 mL/min
Temperature	30°C
Detection	UV at 254 nm

Injection	1.0 μ L (1 mg/mL)
Rt1 & Rt2 (min)	1.88 & 2.18
k1 & k2	2.93 & 3.57
Separation Value	1.22
Resolution Value	2.59