

Reliasil Amino HPLC Column

Chemicals:

Arabinose, mannose, xylose, glucose, galactose

HPLC Condition:

Column: Reliasil Amino HPLC column 4.6x250 mm, 3 μ m.

Part number:


Mobile phase: Acetonitrile: water 75:25 (v/v)

Flow rate: 1 mL/min

Injection volume: 20 μ L

Column temp: room temperature

RI detector

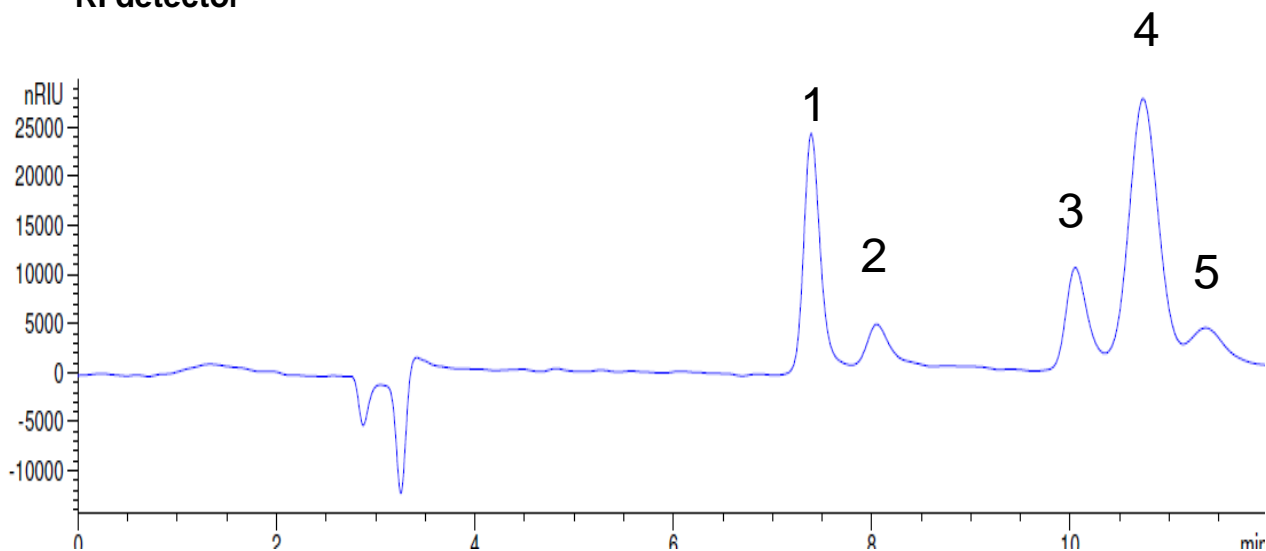

1. Xylose; 2. Arabinose; 3. Mannose; 4. Glucose; 5. Galactose