


Separation of 1-Acenaphthenol on Epitomize CSP-1C by UHPLC

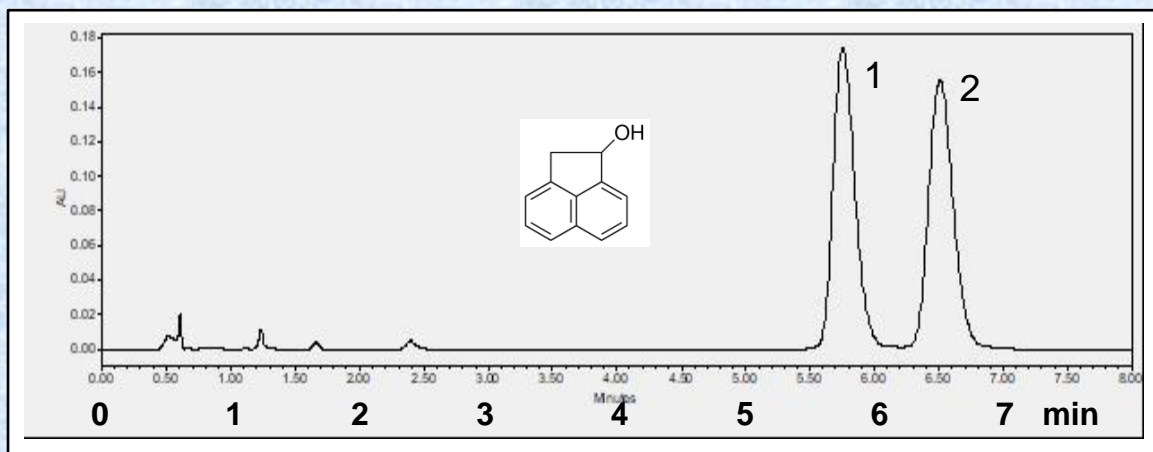

Sample	1-Acenaphthenol [D18724]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1C, 3 micron, 3.0 x 50 mm		
Mobile Phase	97/3 Heptane/2-propanol		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	5.75 & 6.51
Temperature	25 °C	k1 & k2	10.28 & 11.76
Detection	UV at 230 nm	Separation Value	1.14
Injection	0.20 µL (1 mg/mL)	Resolution	2.32