


Separation of Benzoin on Epitomize CSP-1C by UHPLC

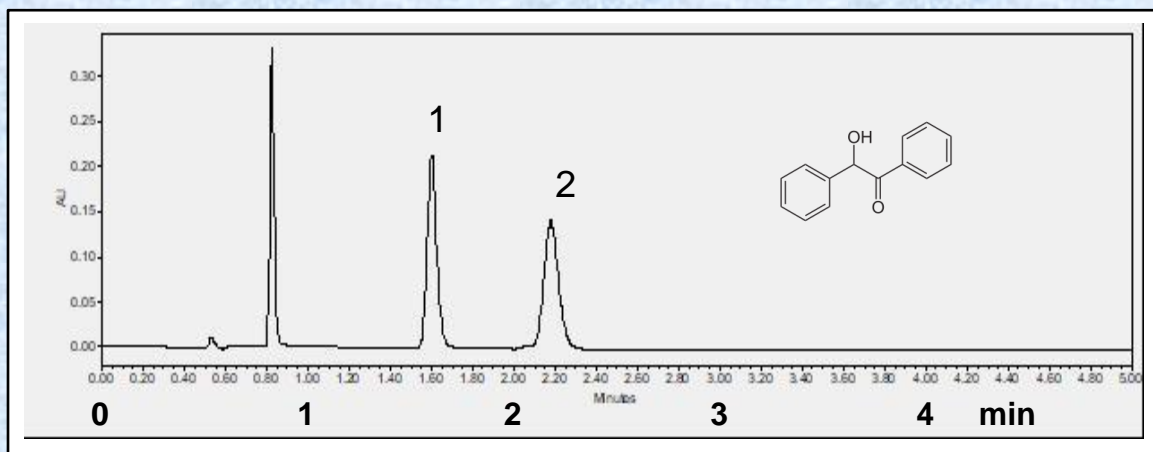

Sample	Benzoin [D24703]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1C, 3 micron, 3.0 x 50 mm		
Mobile Phase	90/10 Heptane/2-propanol		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	1.60 (S) & 2.18 (R)
Temperature	25 °C	k1 & k2	2.14 & 3.27
Detection	UV at 230 nm	Separation Value	1.53
Injection	0.50 µL (1 mg/mL)	Resolution	5.41