


Separation of 1,1'-Bi-2-naphthol on Epitomize CSP-1Z by UHPLC

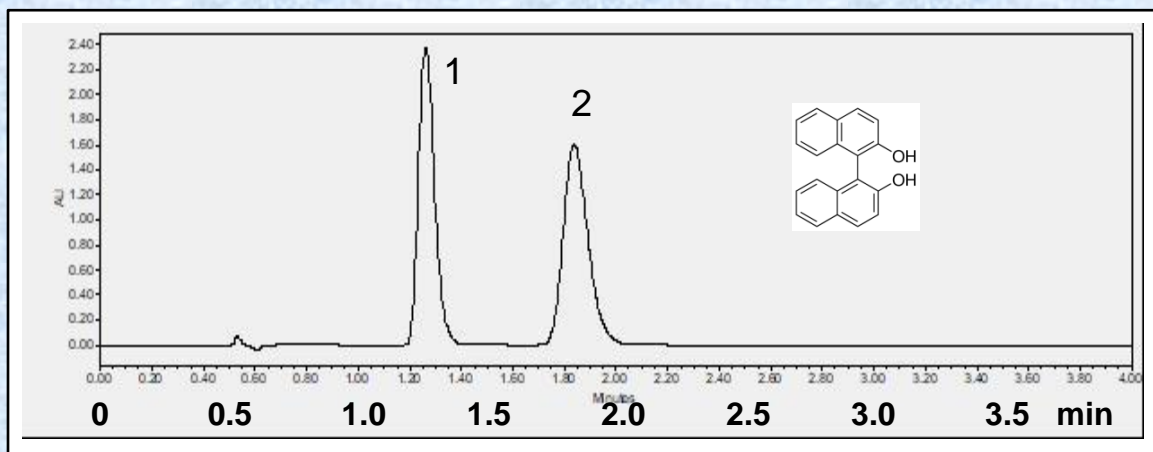

Sample	1,1'-Bi-2-naphthol [D24808]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1Z, 3 micron, 3.0 x 50 mm		
Mobile Phase	90/10/0.1 Heptane/2-propanol/TFA		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	1.26 (S) & 1.84 (R)
Temperature	25 °C	k1 & k2	1.44 & 2.56
Detection	UV at 230 nm	Separation Value	1.77
Injection	0.50 µL (1 mg/mL)	Resolution	4.01