


Separation of Flocinopiperidol on Epitomize CSP-1K by UHPLC

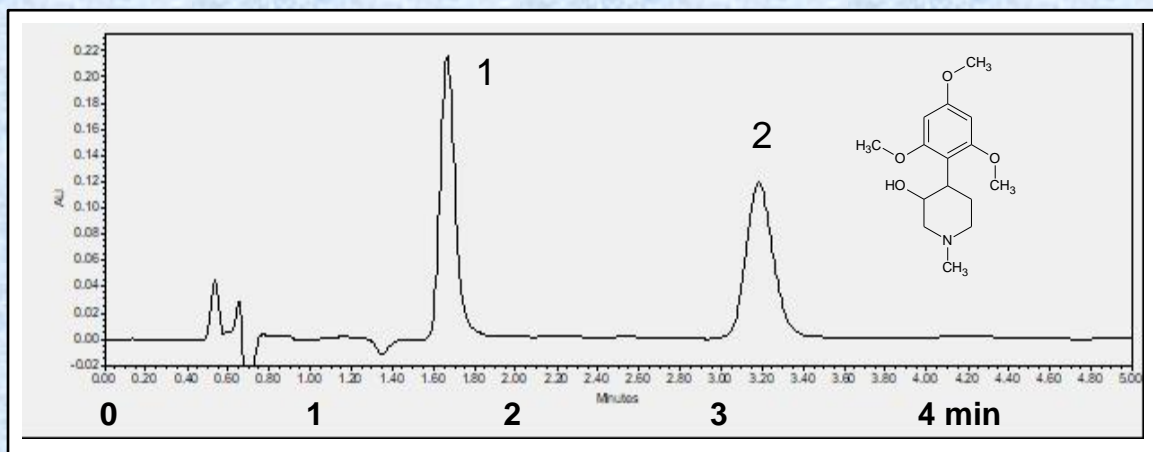

Sample	Flocinopiperidol [D26383]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1K, 3 micron, 3.0 x 50 mm		
Mobile Phase	80/20/0.1 Heptane/2-propanol/DEA		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	1.67 & 3.18
Temperature	25 °C	k1 & k2	2.25 & 5.21
Detection	UV at 230 nm	Separation Value	2.32
Injection	0.50 µL (1 mg/mL)	Resolution	8.35