


Separation of Ketoprofen on Epitomize CSP-1K by UHPLC

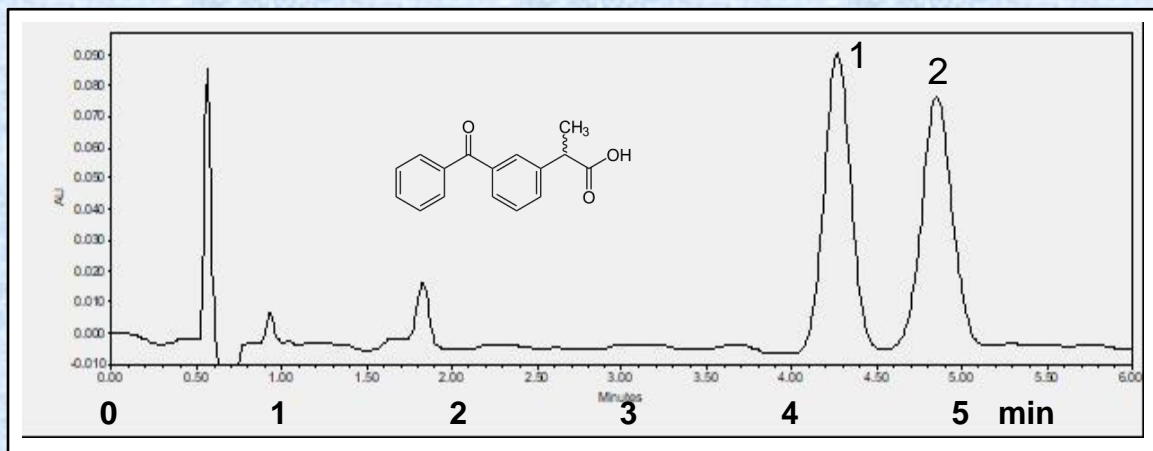

Sample	Ketoprofen [D26169]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1K, 3 micron, 3.0 x 50 mm		
Mobile Phase	95/5/0.1 Heptane/2-propanol/TFA		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	4.27 & 4.85
Temperature	25 °C	k1 & k2	7.32 & 8.46
Detection	UV at 230 nm	Separation Value	1.16
Injection	0.50 µL (1 mg/mL)	Resolution	1.90