


Separation of Miconazole on Epitomize CSP-1Z by UHPLC

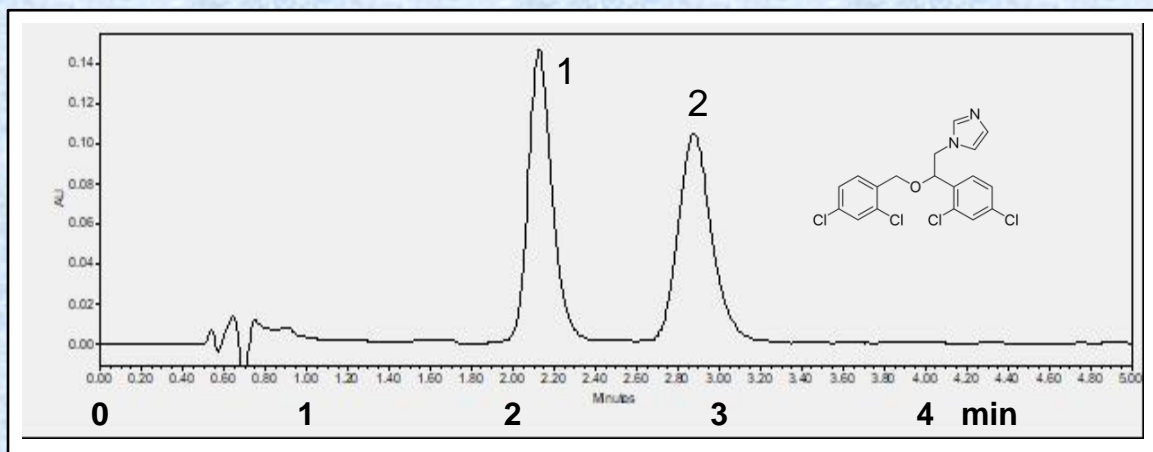

Sample	Miconazole [D23172]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1Z, 3 micron, 3.0 x 50 mm		
Mobile Phase	80/20/0.1 Heptane/2-propanol/DEA		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	2.13 & 2.88
Temperature	25 °C	k1 & k2	3.12 & 4.57
Detection	UV at 230 nm	Separation Value	1.46
Injection	0.50 µL (1 mg/mL)	Resolution	3.03