


Separation of 1-(1-Naphthyl)ethanol on Epitomize CSP-1Z by UHPLC

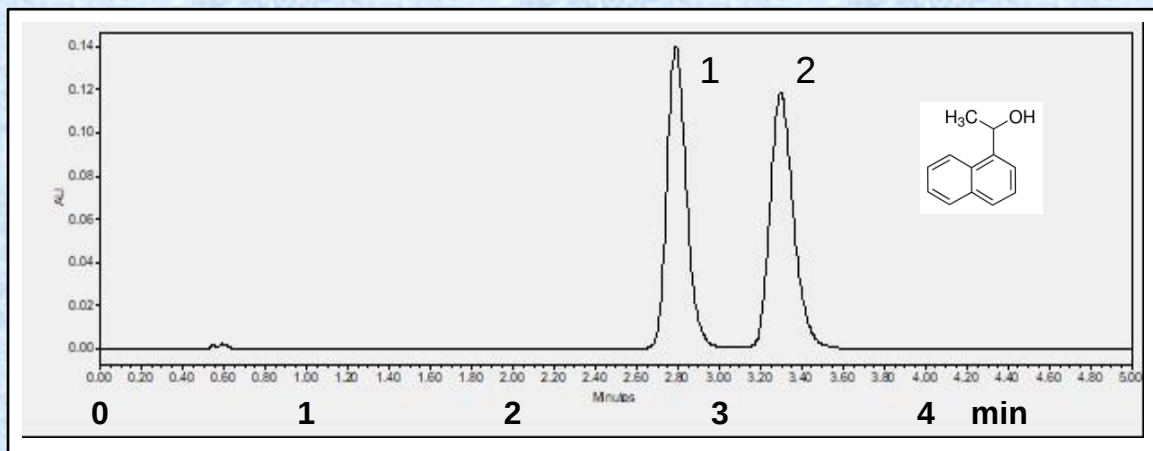

Sample	1-(1-Naphthyl)ethanol [D18939]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1Z, 3 micron, 3.0 x 50 mm		
Mobile Phase	97/3 Heptane/2-propanol		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	2.79 & 3.30
Temperature	25 °C	k1 & k2	4.47 & 5.47
Detection	UV at 230 nm	Separation Value	1.22
Injection	0.20 µL (1 mg/mL)	Resolution	2.73