


Separation of *trans*-Stilbene Oxide on Epitomize CSP-1A by UHPLC

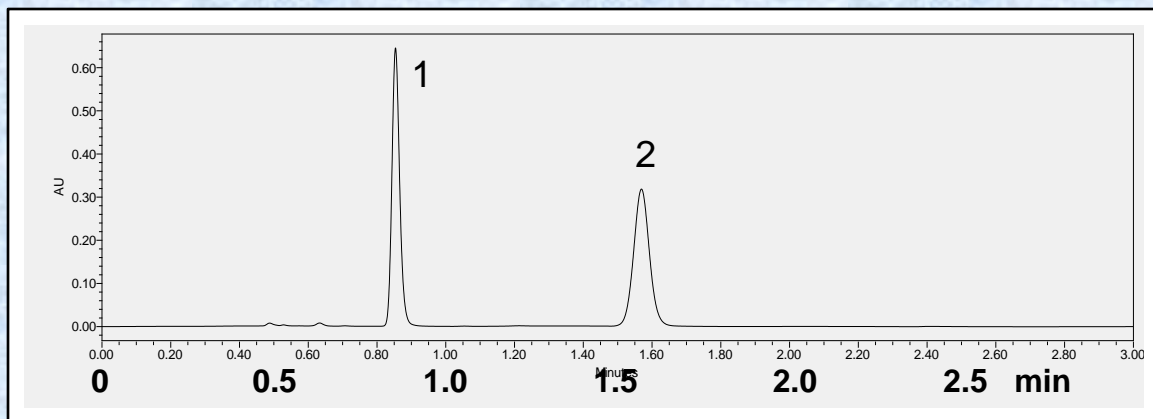

Sample	<i>trans</i> -Stilbene Oxide [D18353]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1A, 3 micron, 3.0 x 50 mm		
Mobile Phase	90/10 Heptane/2-propanol		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	0.85 & 1.57
Temperature	25 °C	k1 & k2	0.80 & 2.31
Detection	UV at 230 nm	Separation Value	2.88
Injection	0.20 µL (1 mg/mL)	Resolution	11.65