


Separation of *trans*-Stilbene Oxide on Epitomize CSP-1C by UHPLC

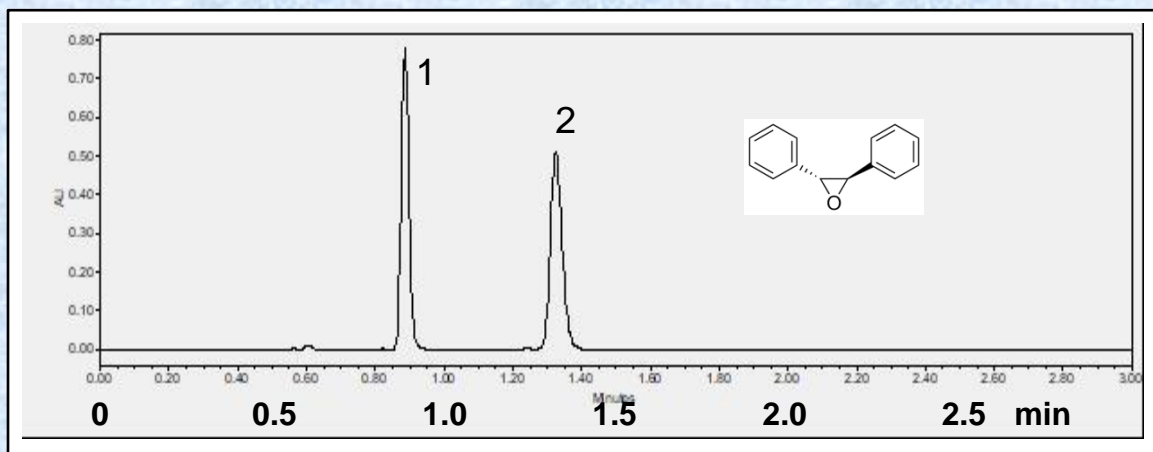

Sample	<i>trans</i> -Stilbene Oxide [D18621]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1C, 3 micron, 3.0 x 50 mm		
Mobile Phase	90/10 Heptane/2-propanol		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	0.89 & 1.33
Temperature	25 °C	k1 & k2	0.74 & 1.60
Detection	UV at 230 nm	Separation Value	2.17
Injection	0.20 µL (1 mg/mL)	Resolution	8.50