


Separation of Tetramisole on Epitomize CSP-1A by UHPLC

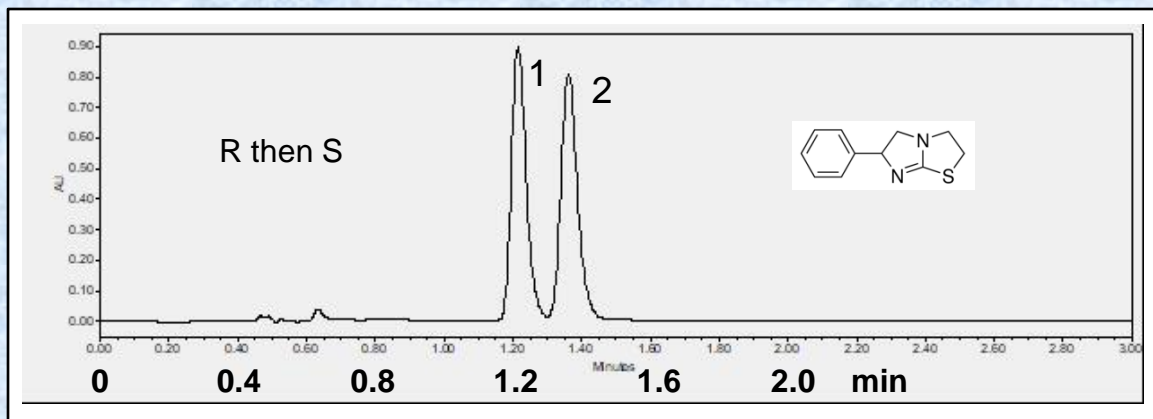

Sample	Tetramisole [D22004]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1A, 3 micron, 3.0 x 50 mm		
Mobile Phase	80/20/0.1 Heptane/2-propanol/diethylamine		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	1.21 & 1.36
Temperature	25 °C	k1 & k2	1.56 & 1.87
Detection	UV at 230 nm	Separation Value	1.20
Injection	1.0 µL (1 mg/mL)	Resolution	1.79