


Separation of Verapamil on Epitomize CSP-1K by UHPLC

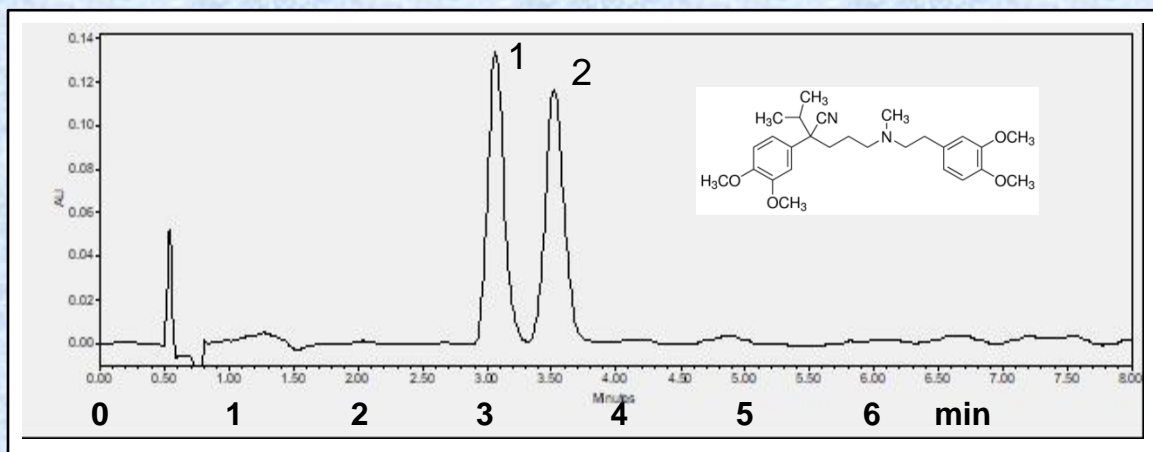

Sample	Verapamil [D26714]		
System	Waters Acquity, H-Class		
Column	Epitomize CSP-1K, 3 micron, 3.0 x 50 mm		
Mobile Phase	90/10/0.1 Heptane/2-propanol/DEA		
Flow Rate	0.50 mL/min	Rt1 & Rt2 (min)	3.06 & 3.52
Temperature	25 °C	k1 & k2	4.97 & 5.86
Detection	UV at 230 nm	Separation Value	1.18
Injection	0.50 µL (1 mg/mL)	Resolution	1.83