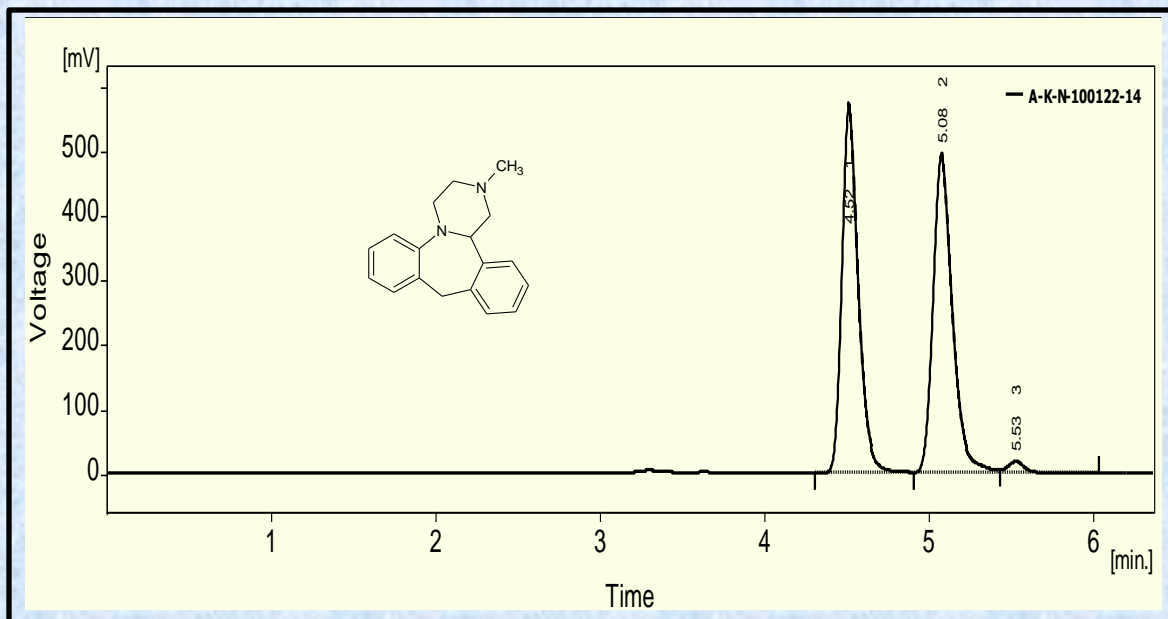


**Mianserin on Epitomize CSP-1K, 5  $\mu$ m**


<b>Column</b>	CSP-1K, 5 micron
<b>Column Size</b>	4.6 x 250 mm
<b>Mobile Phase</b>	98/2 heptane/2-propanol
<b>Flow Rate (mL/min)</b>	1.0
<b>Temperature (°C)</b>	25
<b>Detection</b>	UV at 254 nm
<b>Rt1 &amp; Rt2 (min)</b>	4.52 & 5.08
<b>k1 &amp; k2</b>	0.50 & 0.68
<b>Separation Value</b>	1.38