


Propranolol on Epitomize CSP-1A, 3 μm , SFC-UHPLC

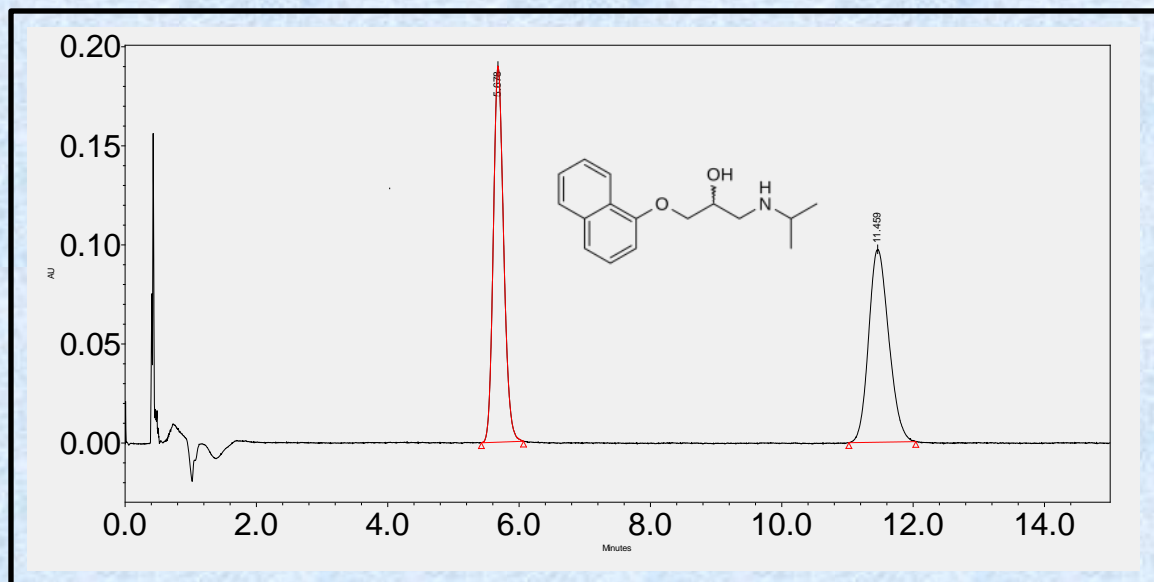

System	Waters Acquity UPC ²
Column	CSP-1A, 3 micron
Column Size	3.0 x 50 mm
Mobile Phase	90/10 CO ₂ /Modifier
Modifier	50/50/0.2 MeOH/IPA/DEA
Flow Rate	0.80 mL/min
Temperature	35°C
Detection	UV at 220 nm

Injection	0.50 μL (1 mg/mL)
Rt1 & Rt2 (min)	5.68 & 11.46
k1 & k2	10.88 & 22.97
Separation Value	2.11
Resolution Value	13.72